EN

Thank you for purchasing our scale with speaker function. This manual will provide instructions for effective and safe use of the product. Please read it carefully before use and store the manual in a safe place for future reference. If you ever pass the product on to others, this manual should be included.

1. Technical data

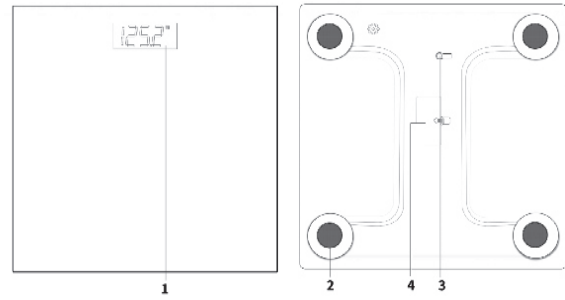
| | |
|-----------------|---------------------------------|
| Display | LCD |
| Product size | 28 x 28 x 2,5 cm |
| Weight capacity | 2,5-180 kg |
| Weight accuracy | 0,1 kg |
| Power | 2x AAA batteries (not included) |

2. Safety instructions

- Do not immerse in water.
- Make sure your feet are dry before stepping on the scale.
- Do not stand on the edge of the scale or jump on it.
- Do not drop the device or drop objects on it. This may damage the sensors.
- Do not overload the scale (max 180 kg).
- Place the device on a hard, level and dry surface.
- Do not attempt to repair, disassemble or alter the product yourself.

3. Description of parts

1. LCD display
2. Non-slip feet
3. Control button
4. Battery cover

**4. Operation****Inserting batteries**

1. Remove the battery cover on the back of the scale.
2. Insert 2x AAA batteries, with regards to polarity (+ and -).
3. Replace the battery cover.

Note: Incorrect polarity may damage electronic components.

Using the body scale

1. Place the scale on a hard, level and dry surface.
2. Step briefly onto the scale with dry feet to activate it.
- Wait until the display shows 0,0 kg.
- Stand still for an accurate measurement.
3. The scale will announce and display your weight.
4. It switches off automatically after about 10 seconds.

Note:

- For the most accurate results, weigh yourself at the same time every day.
- Also wake up the scale to recalibrate it if it has been moved.

Control button

Wake up the scale and press the control button, located on the back of the scale to:

• Change languages:

English Dutch French German

En NL Fr DE

• Enable/disable audio: press and hold for 3 seconds

Voice and sounds enabled

on

Voice and sounds disabled

off

Error indication

When an error has occurred, the scale will give an audio signal and display the following symbols:

| | |
|----|--|
| LO | The battery power is running low. Please replace the batteries. |
| OL | The scale is overloaded. Please step off to avoid possible damage to the device. |

5. Disposal

Dispose of the article and the packaging materials in accordance with current local regulations. This product is subject to the provisions of the European Directive 2012/19/EU.

The crossed out wheeled bin symbol means that this product shall not be disposed of with normal household waste. Electronic and Electrical Equipment not included in the selective sorting process are potentially dangerous for the environment and human health due to the presence of hazardous substances. Please dispose of responsibly at an approved waste or recycling facility.

The CE marking is a certification mark that indicates conformity with health, safety and environmental protection standards for products sold within the European Union. The CE marking is the manufacturer's declaration that the product meets the requirements of the applicable EC directives.

Recycling helps to reduce the consumption of resources and reduce the environmental impact. If you wish to discard the packaging, comply to the environmental legislation in your country.

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