 SILVERGEAR

**EYE MASSAGER**

RELAX AND RECOVER

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Thank you for purchasing the Silvergear Eye massager. This manual will provide instructions for effective and safe use of the device. Please read it carefully before use and store the manual in a safe place for future reference.

### 1. Package content

- Eye massager
- USB cable
- Manual

### 2. Technical data

Adapter: AC 100-240V DC 5V

Battery: 1200mAh

Weight: 375g

Belt Length: 38cm

Belt Width: 3.2cm

### 3. Safety instructions

- For intended use only, as described in this manual.
- Do not use while driving, walking or operating appliances.
- Do not place, store or use the product near water.
- Do not use the product in a hot or humid environment.
- Before use, remove contact lenses.
- Stop use immediately in case of irritation or pain.
- Not recommended for children, pregnant women, or patients with eye disease.
- Prior to use, consult your general practitioner in case of heart conditions, vascular disorders or high or low blood pressure.
- When undergoing treatment for any eye condition, consult your doctor prior to use.
- Do not attempt to repair, disassemble or alter the product yourself.

### 4. Description of parts

1. Power button

To turn the device ON / OFF and to switch between the different modes.

2. ⏮ -


Short press: play the previous song - Long press: turn the volume DOWN

3. ⏭ +

Short press: play the next song - Long press: turn the volume UP

4. Indicator light

5. Micro USB charging port
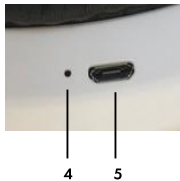
### Display:

- Bluetooth
- Music (plays/pause)
- Air pressure
- Vibration
- Hot compress
- Battery level

### 5. Operation Charging

Before first use, make sure you charge the device fully for 3-5 hours. Connect the Eye massager with the supplied USB cable to a power supply (computer / power bank / mains). The charging indicator light will glow red during charging. The Eye massager can be used for 8-10 sessions on a single charge.

### Using the Eye massager

First, remove contact lenses for safety and any earrings, hair pins etc. for more comfort.

1. Before use, unplug the device.
2. Press and hold the power button for 3 seconds to turn on the device. Bluetooth is automatically activated. Search on your mobile phone for the Eye massager (model: Eye Massager 2465) and make a connection.
3. Repeatedly press the power button to switch between these 5 modes:
  - Integrated mode (air pressure + vibration + hot compress + music)
  - Clear mode (air pressure + hot compress + music)
  - Sleeping mode (air pressure)
  - Vitality mode (air pressure + vibration + music)
  - Comfort mode (hot compress + music)
4. Press and hold the power button for 3 seconds to turn off the device. It also turns off automatically after 15 minutes.

### 6. Maintenance

Clean the device with a slightly damp towel and wipe it dry with a soft cloth. Store it in a dry, cool place away from direct sunlight.

**Note:** Do not use scouring pads or an abrasive detergent to clean the device.

### 7. Disposal

The crossed out wheeled bin symbol on this product is to remind you that electrical products should not be disposed of with household waste. Please dispose of it responsibly at an approved waste or recycling facility. Check with your local authority for recycling advice.

Dispose of all packaging materials in an environmentally friendly manner.

CE marking is a certification mark that indicates conformity with health, safety and environmental protection standards for products sold within the European Union. The CE marking is the manufacturer's declaration that the product meets the requirements of the applicable EC directives.

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