Microscope


Science experiment

Scientific Microscope Set

EN

MICROSCOPE PARTS

1. Ocular
2. Tube
3. Revolving turret (with 3 lenses)
4. Objective lenses
5. Stage
6. 4-colour filter
7. Light unit
8. Knob
9. Arm
10. Specimen clamp
11. ON/OFF switch
12. Base

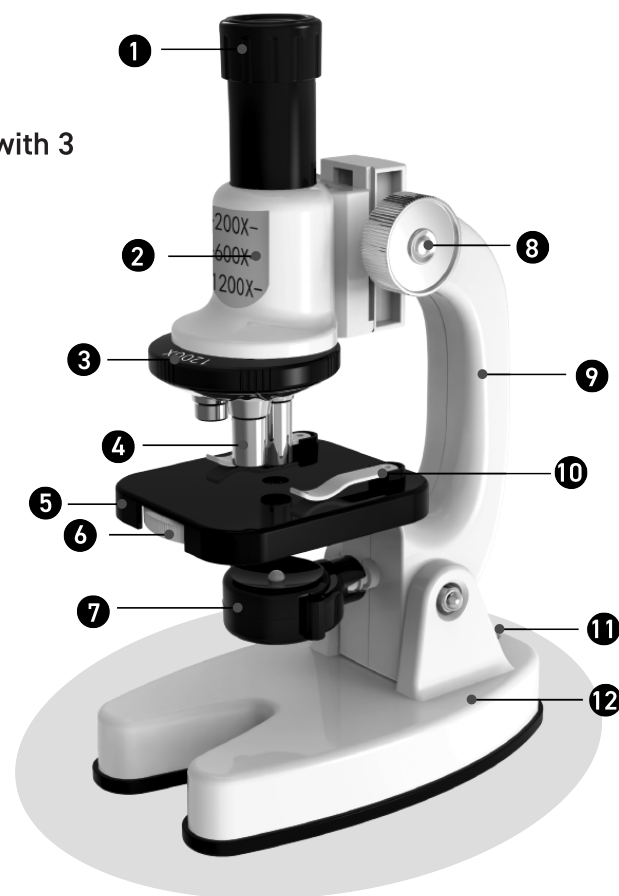
HOW TO USE THE MICROSCOPE

1 Turn on the LED light in the light unit with the LED switch. Make sure to add 2x 1.5V "AA" batteries to the base of the microscope.

2 Alternatively, use the mirror on the back of the lighting unit (no batteries required). Tilt and rotate the unit so that external light is reflected through the round hole in the centre of the stage.

4 Turn the knob until the lenses are at their highest position. Begin observation with the lowest magnifying objective and work up to the 1200x objective. Remember, the longer the objective, the greater the magnification.

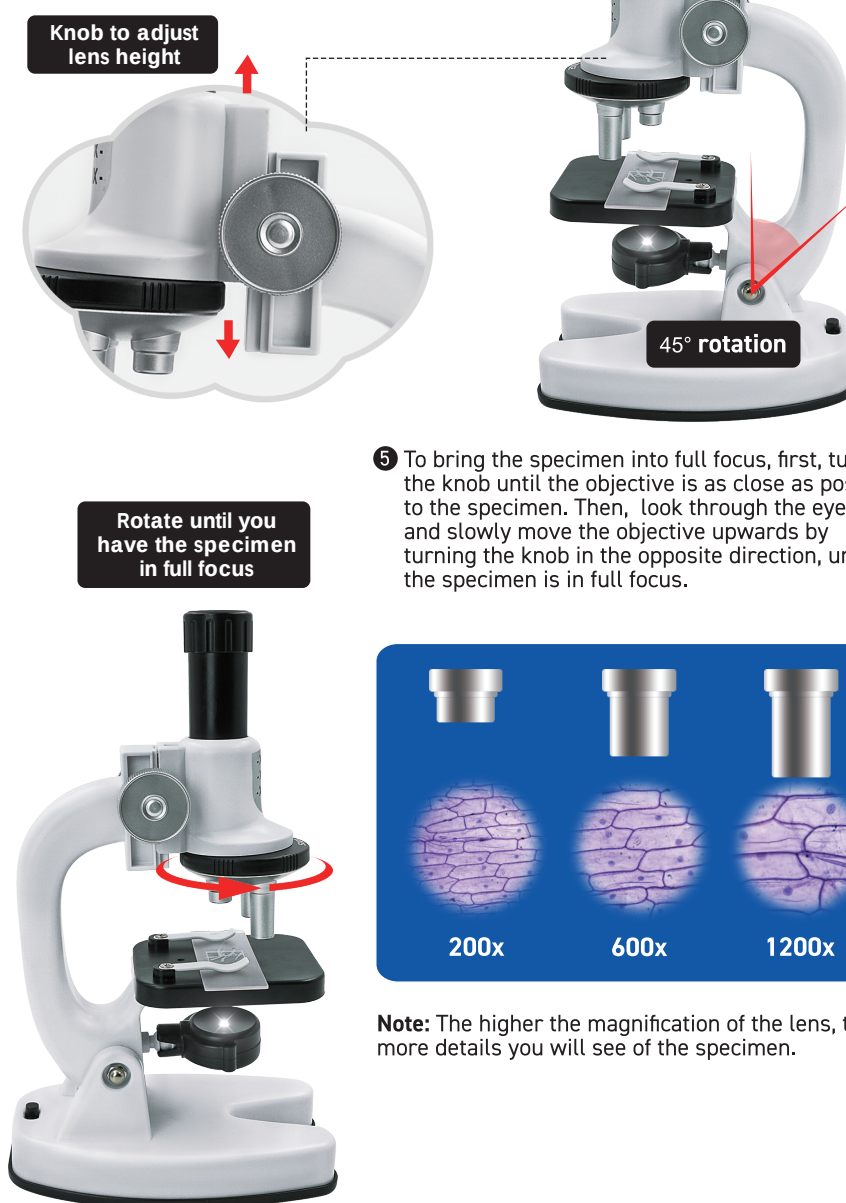6 The 4-colour filter has been installed at the bottom of the stage. rotate the filter to switch between colours and experiment with the different observation effects.

Green Blue Yellow Red

INSTALLATION OF PHONE HOLDER

1 Remove the protective film on one side of the transparent stickers and carefully stick them to the phone holder. Place the phone holder on top of the eyepiece and remove the remaining protective film.

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